

WIDE $PIXL_{2(1)\times15-MPX3}^{@}$ Datasheet

Model No.: WxCM3x-Xxx210721





General features



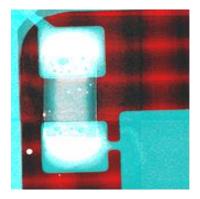


Illustration of multichannel "color" radiographs where different materials are identified and imaged in different colors

The large area imaging detector **WIDE PIX** $L_{2(1)x15-MPX3}$ with resolution of 512 (256) x 3840 pixels is composed of Medipix3 hybrid detector electronics tiles. Each tile (256 x 256 pixels) is attached to a silicon or CdTe sensor. Therefore, the whole area of the **WIDE PIX** $L_{2(1)x15-MPX3}$ device is fully sensitive and there are no gaps between sensor tiles. Each pixel has two integrated 12-bit digital counters and two energy discrimination thresholds. The counters store number of registered particles, e.g. X-ray photons, with energy above the appropriate threshold. Both counters can be joined to a single 24-bit counter providing enhanced dynamic range. The particle counting principle eliminates any additional noise generated by the sensor or electronic readout. It allows acquiring X-ray images with very high contrast and wide dynamic range. Therefore, even low contrast structures such as plastic or soft tissue are easily detectable in X-ray images.

Both devices are suitable for CT scanners, which can take advantage of large sensitive area without any gaps. The **WIDE**PIX L_{1x15} - MPX3 variant moreover supports a hardware-based Time-Delayed-Integration mode for online (continuous) scanning applications.

The energy discrimination thresholds of Medipix3 technology allow spectral X-ray imaging. Different materials in an inspected sample could be then identified based on their spectral X-ray attenuation properties. Energy spectra could be measured typically from 5 keV upwards.

The Charge Summing Mode implemented in the pixel electronics provides hardware-based correction of signal cross talk between pixels. This further considerably improves the detector spectral response and therefore also quality of spectra measured in individual pixels.

The camera is connected to a computer via three ethernet cables.

Main Features

•	Readout chip type	. Medipix3
•	Pixel size	. 55 x 55 μm*
•	Sensor resolution	. 512 (256) x 3840 pixels
•	Dynamic range in one frame	. 12-bit / 24-bit ¹
•	Dark current	. none
•	Interface	.3x RJ45 1Gbit/s ethernet
•	Maximum frame rate	. up to 80 (170) fps¹
•	Dimensions	. 280 x 153 x 42 mm
•	Weight	3500 g



¹ Depends on operation mode.

^{* 55} x 110 μ m at the edges and 110 x 110 μ m at the corners



Device parameters

Operating conditions

Symbol	Parameter	Min	Тур	Max	Units	Comment
T _A	Temperature Range	0	30	40	°C	
Φ	Humidity		55	60	%	Not condensing
	IP Code		IP50			

Family parameters

 $T_A = 25$ °C, $V_{CC} = 24V$

Symbol	Parameter	WidePIX L 1x15 - MPX3	WidePIX L 2x15 - MPX3	Units	Comment
Vcc	Supply Voltage	20/24/26		V	Min/Typ/Max
Icc	Supply Current (V _{CC} = 24V)	1000/2000	2400/4800	mA	Typ/Max
P	Power dissipation	27/54	48/96	W	Typ/Max
Α	Sensor Area	214.5 x 14.1	214.5 x 28.2	mm	
	Detector Resolution	3840 x 256	3840 x 512	Pixels	
f	Frame Rate ¹	50	20	fps	
T _{READ}	Readout Time ²	20	50	ms	
m	Weight	3300	3500	g	

Modes of readout chip operation

Туре	Mode	Bit depth	Description
	SPM-1CH	12/24 bit/frame	Single Pixel Mode using one counter: Every pixel works independently of its neighbors. One energy threshold (energy channel) is available. 1 output image: Number of events per pixel
Frame (reading	SPM-2CH	12 bit/frame	Single Pixel Mode using both counters: Every pixel works independently of its neighbors. Two energy thresholds (energy channels) are available. 2 output images: Number of events per pixel
all pixels)	CSM	12/24 bit/frame	Charge Summing Mode: The charge from 4 adjacent pixels is summed and is assigned to the pixel with the largest charge deposition. The event is counted only if the sum of signals exceeds the energy threshold. 1 output image: Number of events per pixel

All modes can be operated at three ranges: Broad / Narrow / Super Narrow³



¹ Operating parameters: Shutter time=1ms, Mode = CSM or SPM-1Ch 12bit resolution.

² During Readout time (or Dead time), no charge is collected from the sensor.

³ Except CSM mode Silicon sensor.



Energy range and resolution

Typical values for 300 μm Silicon sensor, $T_A = 20$ °C.

Range	Mode	Min Energy [keV]	Max Energy [keV]
Super Narrow	SPM	2.5 - 4.0	100 - 110
Narrou	SPM	2.5 – 4.0	160 - 180
Narrow	CSM	6.0	-
Broad	SPM	5.0 – 6.5	260 - 290
Broad	CSM	8.0	-

Typical values for 1000 μm CdTe sensor, $T_A = 20$ °C

Range	Mode	Min Energy [keV]	Max Energy [keV]
Cupor Norrow	SPM	5 – 10	150 – 170
Super Narrow	CSM	5 – 10	150 – 160
Narrow	SPM	5 – 10	170 – 190
INATIOW	CSM	5 – 10	260 – 270
Drood	SPM	5 – 10	220 – 240
Broad	CSM	5 – 10	460 – 470

^{*} resolution in Sigma of gaussian fit.

Sensor parameters

Sensor material


Symbol	Parameter Si		CdTe	Units	Comment	
	Thickness	300	500	1000	μm	
V _{bias}	Bias Voltage	200	300	-450	V	Max
	Typical detectable energy range for X-rays*	up t	o 60	up to 600	keV	See chart below
	Pixel size	55 x	55**		μm²	

^{*} to get true detector response, detectable energy and quantum efficiency of sensor chip has to be combined with energy range of readout chip (see chapter "Energy range and resolution")




^{**} Pixels on tile borders are 2.5 times larger in one direction. The corner tile pixels are 2.5 times larger in both directions.



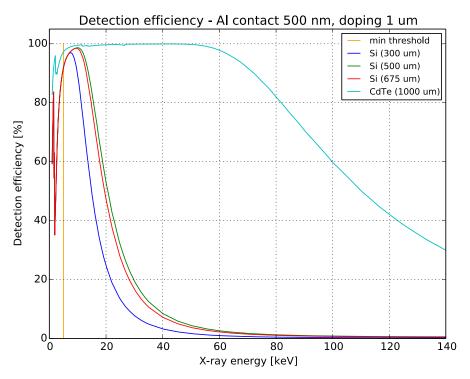

Device description

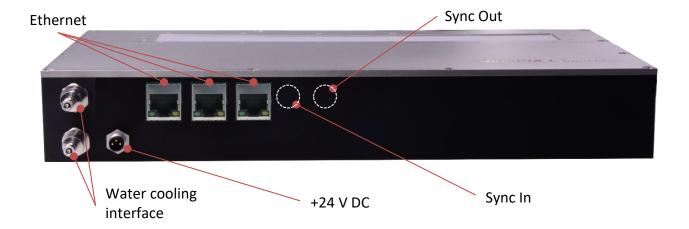

Ethernet connectors

3x RJ45 1Gbit/s ethernet connectors.

+24 V DC connector

Main power supply (via standard M8 connector with 3 female contacts) Connect after plugging USB connector.





Water cooling interface

It is mandatory to cool down detector when in operation. **WIDE** PIX $L_{2(1)x15-MPX3}$ uses water connectors that allow for quick disconnection/reconnection. Mating connector is included as standard accessories and has to be attached to 4x6 mm plastic hose.

Synchronization interface (optional)

Two 4-pin M8 connectors (female for outputs and male for input) serve as synchronization interface, allowing to synchronize **WIDE** PIX $L_{2(1)x15}$ - MPX3 detector with external processes. Four signals are available:

- Ready in measurement is not possible, when signal at logical zero
- Trigger in logical zero starts shutter (measurement)
- Ready out logical one if device is ready to for new shutter
- **Trigger out** mirrors shutter (logical zero when shutter is active)


All signals are TTL compatible and 5V tolerant. For detailed description see **Synchronization Guide**.

-	nc. Outputs 8-4Female)	-	nc. Inputs 8-4Male)
Pin	Pin Signal		Signal
1	Gnd	1	Gnd
2	Trigger Out	2	Trigger In
3	Ready Out	3	Ready In
4	Reserved	4	Reserved






Mechanical dimensions

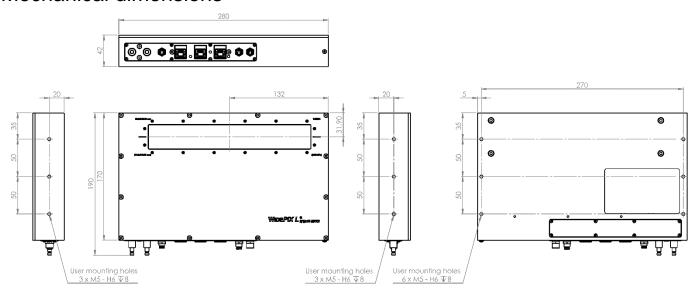

Following drawings are for rear connections

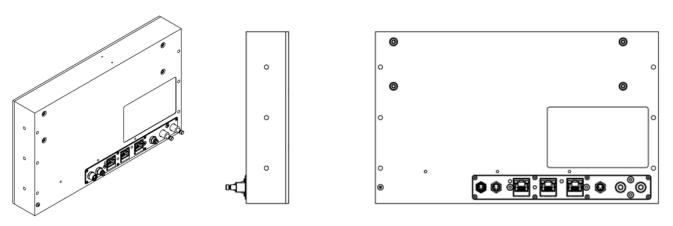

All dimensions are in mm.


* Sensitive surface distance from top of the box may vary depending on actual sensor thickness.






Sensitive area


Sensitive area for models wucm3B-XxxYYMMDD (single upper row)

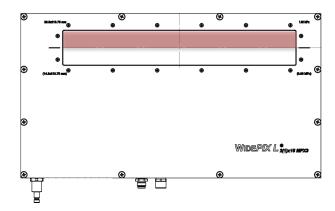

Sensitive area for models wlcm3B-XxxYYMMDD (single lower row)

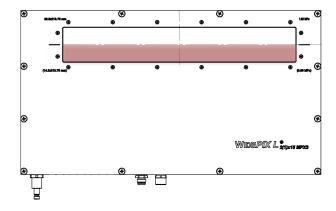


Sensitive area for models w2cm3B-XxxYYMMDD (two rows)

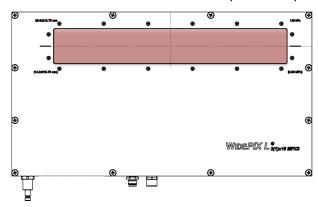





Model Number Codes

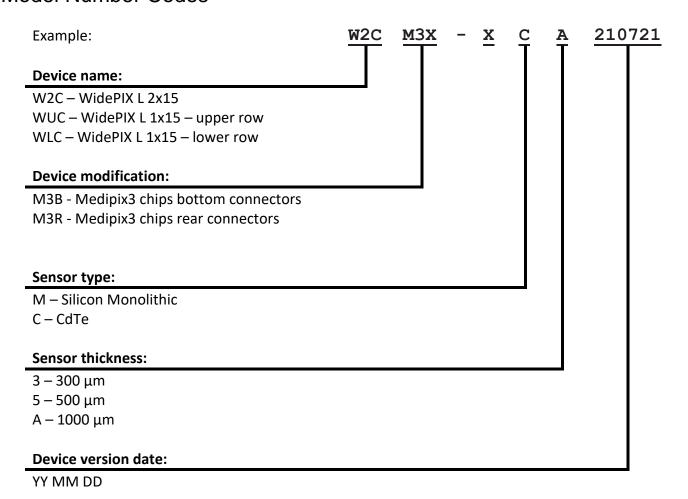

Release history

Date	Changes		
28/07/19 Preliminary version			
26/06/21	ETH version		
23/03/14	New drawings and corrected versions		





Warning

Do not touch sensor surface!

Instructions for safe use

To avoid malfunction or damage to your **WIDE** $PIX L_{2(1)x15-MPX3}$ please obey the following:

- Do not expose to water or moisture **WIDEPIX** $L_{2(1)x15}$ MPX3 is dust protected only.
- Do not open WIDEPIX L_{2(1)x15} MPX3 case. Detector wire-bonding connections may be irreversibly damaged.
- Do not operate detector when not properly water cooled. Otherwise detector temperature may rise above the specified range.

Copyright

ADVACAM s.r.o. Tel: +420 608 605 533

U Pergamenky 1145/12 Email:

170 00 Praha <u>sales@advacam.cz</u>

Czech Republic <u>www.advacam.com</u>

